



Single Use Rotatable Clip Fixing Device

QuickClip 2[™] QuickClip 2 Long[™]

A versatile product line and rotatable mechanism provide
a wider range of hemostatic options than ever before

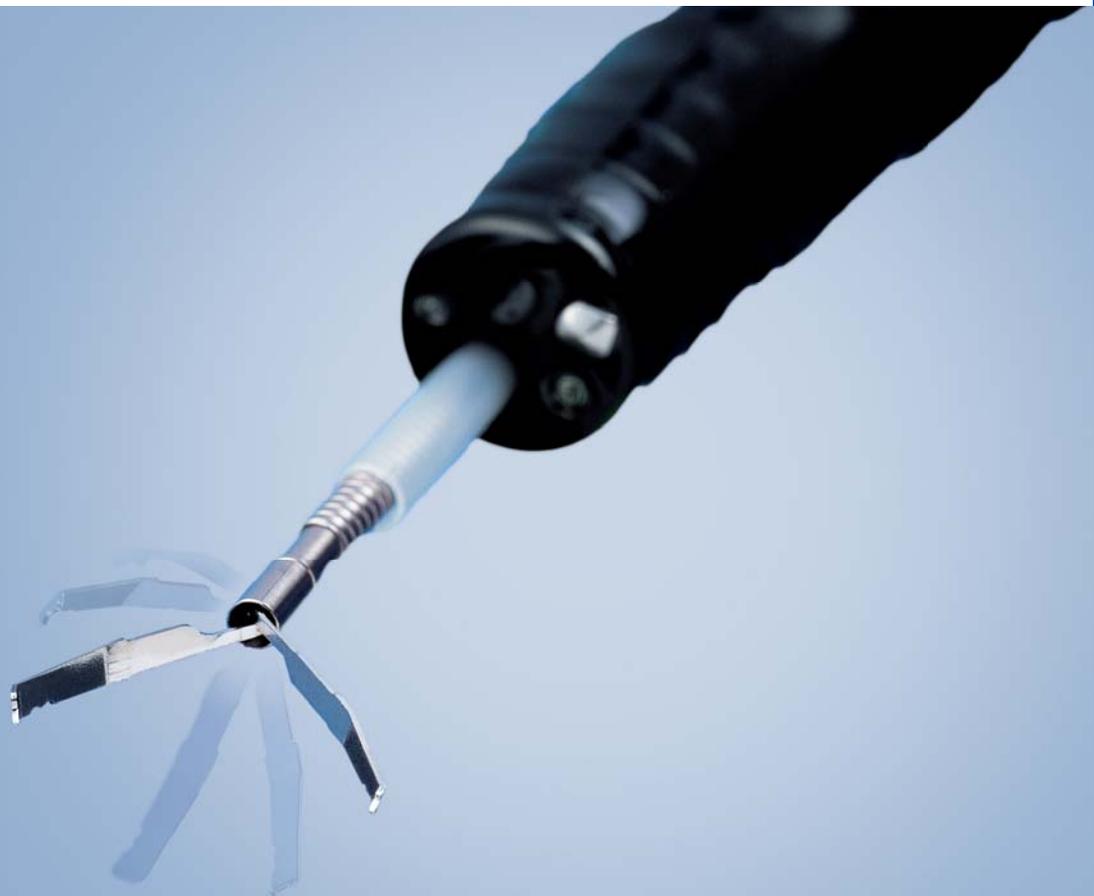


Already featuring models with exceptional rotation abilities, OLYMPUS' versatile line of QuickClip2™ and QuickClip2 Long™ devices now includes new models dedicated to the small intestine

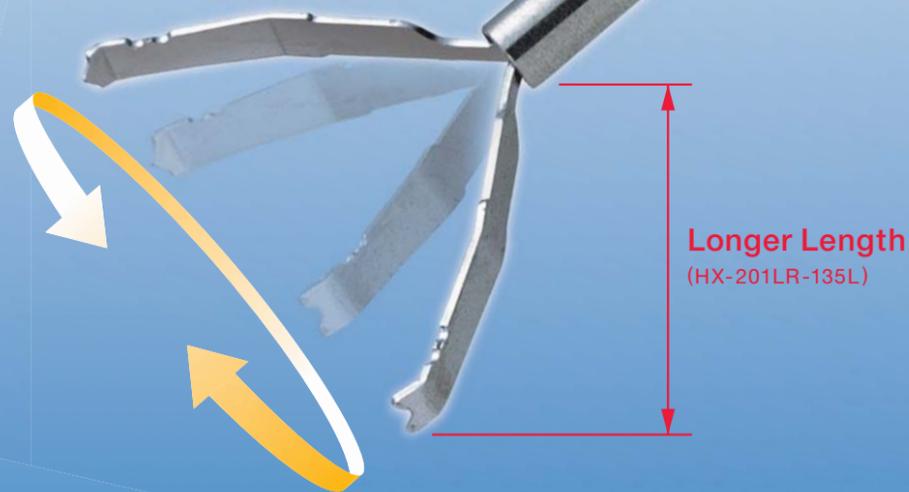

As the creator of endoscopic clipping — widely regarded as one of the most reliable hemostatic techniques in use today, OLYMPUS is committed to maximizing the efficacy of this powerful technique. Our versatile line of QuickClip2™ already includes a proprietary design for smooth and consistent rotation as well as models with various opening widths for differing clinical needs. Now we've added models dedicated to the small intestine, giving you an even wider range of choices and ensuring your hemostatic clipping procedures will be fast and effective no matter what the situation.

QuickClip2™ QuickClip2 Long™

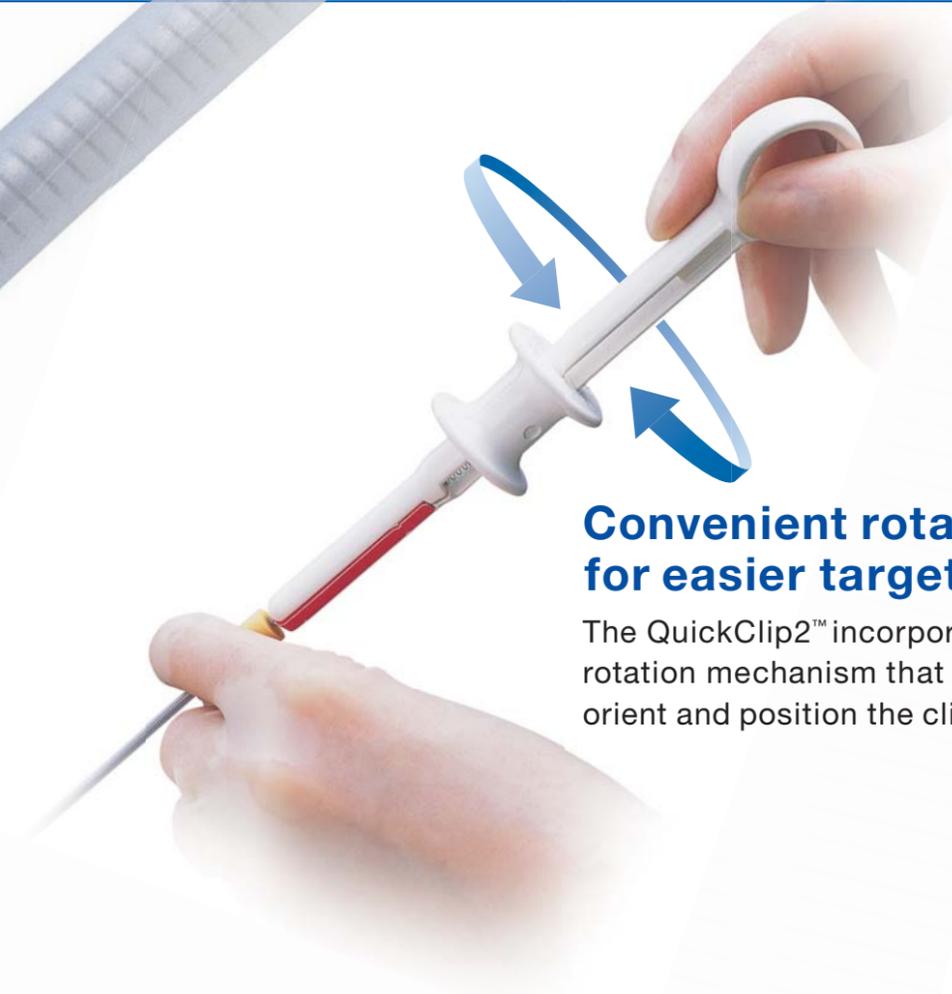


Convenient rotation mechanism for easier targeting

The QuickClip2™ incorporates an exclusive rotation mechanism that enables you to more easily orient and position the clip at the site.


Longer clip length for expanded clearance

As the name suggests, the QuickClip2 Long™ is longer than its predecessor, the original QuickClip2™. The new design gives a much greater opening width and enables more tissue to be grasped at one time — a perfect compliment for the QuickClip2™ family.

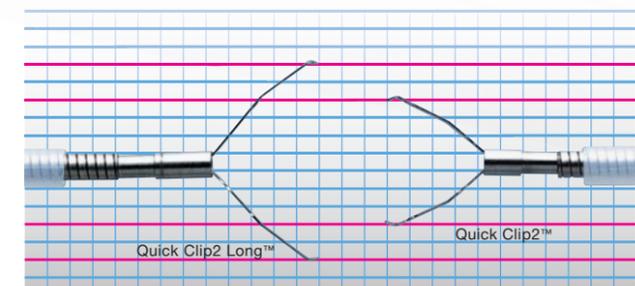

Simple and reliable hemostasis

Endoscopic clipping, considered to be one of the simplest and most reliable hemostatic techniques, mechanically binds blood vessels, minimizing the risk of re-bleeding and damage to the surrounding tissue, while providing instant visual feedback to confirm placement.

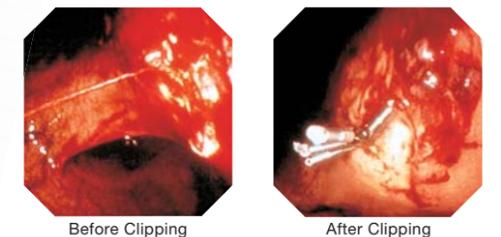

Wider, more versatile lineup

The QuickClip2™ and QuickClip2 Long™ lineup features models with a rotation mechanism designed for use in the upper gastrointestinal tract and models with longer prongs for use in the colon. Now the lineup has been expanded with the addition of new models specifically for use in the small intestine.

Single Use Rotatable Clip Fixing Device

Specifications

Model	Upper GI		Lower GI		Small intestines
	HX-201LR-135 A/B	HX-201LR-135L A/B	HX-201UR-135 A/B	HX-201UR-135L A/B	HX-201YR-135
Working Length	1650 mm		2300 mm		2700 mm
Quantity	A set : 20 pieces B set : 5 pieces		A set : 20 pieces B set : 5 pieces		5 pieces
Compatible Channel	2.8 mm		2.8 mm		2.8 mm
Maximum Diameter	2.6 mm		2.6 mm		2.6 mm
Clip Length	Standard	Long	Standard	Long	Standard
Clip Angle	135°		135°		135°


Intended Use

This instrument has been designed to be used with an OLYMPUS endoscope for endoscopic clip placement within the gastrointestinal (GI) tract for the purpose of


- (1) endoscopic marking,
 - (2) hemostasis for
 - (a) mucosal / sub-mucosal defects < 3 cm,
 - (b) bleeding ulcers,
 - (c) arteries < 2 mm,
 - (d) polyps < 1.5 cm in diameter,
 - (e) diverticula in the colon,
 - (3) as a supplementary method, closure of GI tract luminal perforations < 20 mm that can be treated conservatively.
- Do not use these instruments for any purpose other than their intended use.

QuickClip2™/QuickClip2 Long™: Operating procedure


1 Once the device is seen at the distal end of the endoscope, remove the red stopper. The red stopper prevents the clip from extending out of the plastic sheath during insertion. Pull the yellow tube joint towards the handle to extend the clip from the plastic sheath.

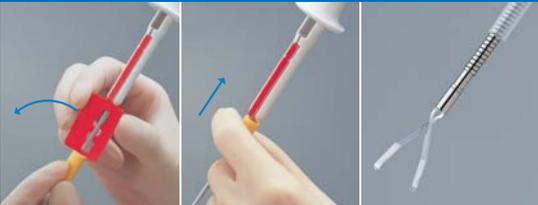

2 Before approaching the targeted bleeding site, it is important to open the clip to its maximum width. To do this, gently pull the slider a very short distance towards the thumb ring.

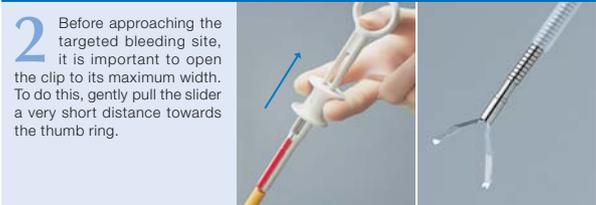

3 While holding the yellow tube joint, turn the handle slowly until the clip is in the proper position for clipping the target site.

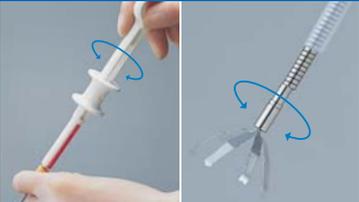

4 While pressing the clip against the target site, pull the slider towards the thumb ring to grasp the target site and deploy the closed clip. Remove the device from the endoscope.



Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.



OLYMPUS MEDICAL SYSTEMS CORP.
Shinjuku Monolith, 2-3-1 Nishi-Shinjuku, Shinjuku-ku, Tokyo 163-0914, Japan

For a complete listing of sales and distribution locations visit:
www.olympus.com