

OLYMPUS

Your Vision, Our Future

ENDOUROLOGY

High-Flow Nephroscope

OES Pro



GO WITH THE PERFECT FLOW

INTRODUCING THE OES PRO

By combining the latest hybrid-lens technology with a high-flow sheath design, the new Olympus percutaneous nephroscope offers unsurpassed

clear image quality. In addition, the extra long working length is not only perfectly suited for obese patients but can be used as a cystoscope as well.

- High flow sheath design
- 30° direction of view
- Crystal clear image quality due to hybrid-lens technology
- 70% more flow with 4 mm instruments compared to predecessor model
- 230 mm working length




The new three-nail grasper WA63817A provides a lumen for laser treatment while holding the stone in a secure position.




Thanks to a special designed obturator, the nephroscope can be used as a compact cystoscope as well.

NEPHROSCOPE

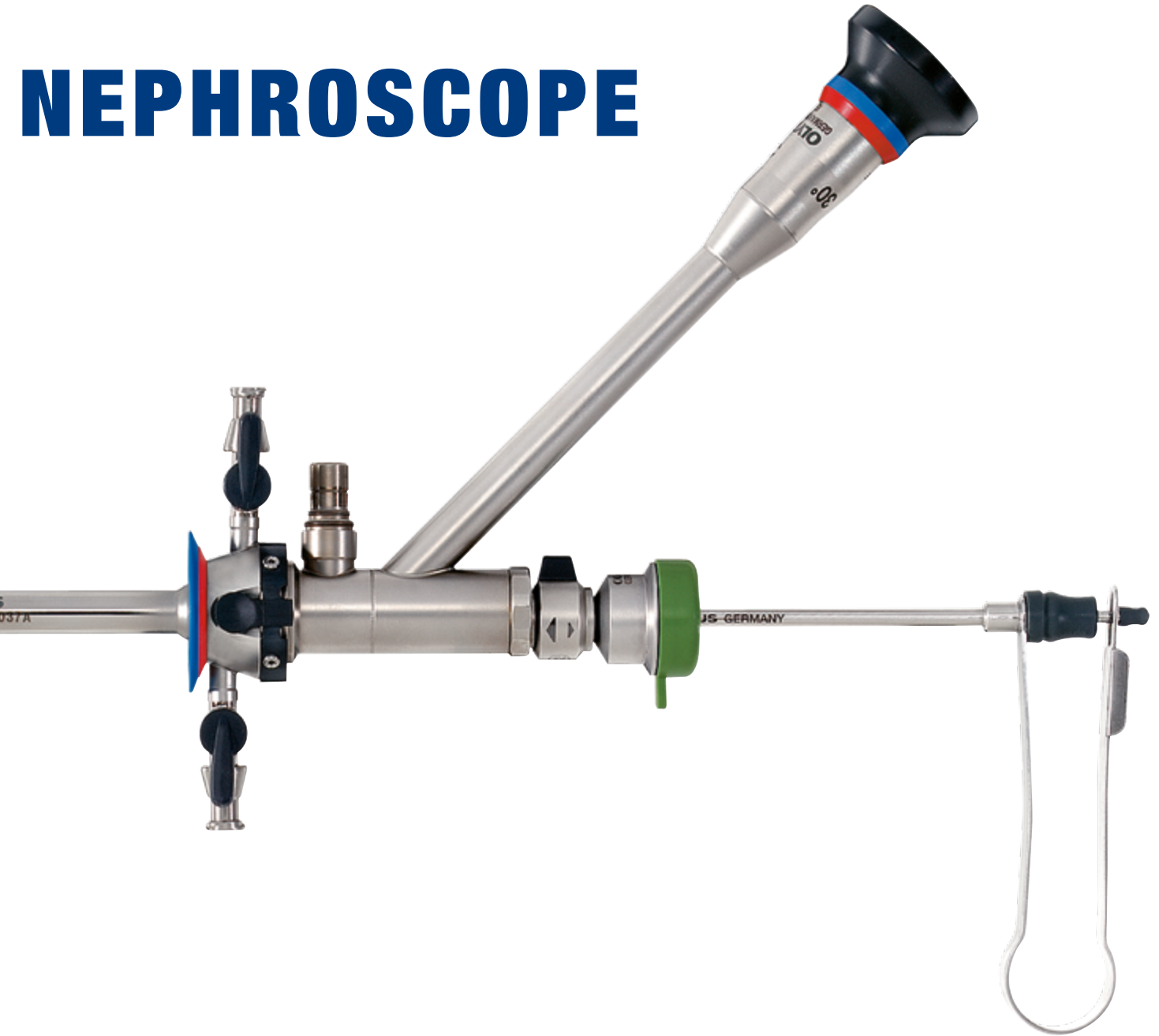


nephroscope

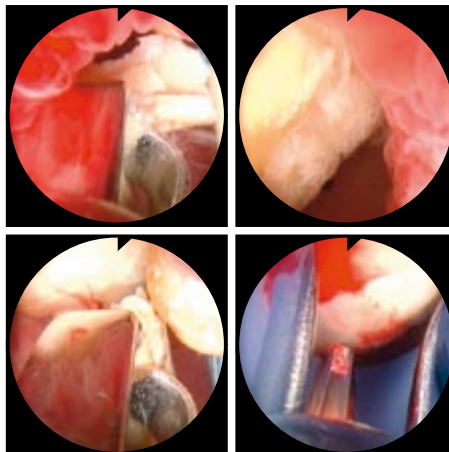

The optical system with 30° direction of view allows early viewing of instrumentation when being advanced through the working channel. In addition, the hybrid-lens technology provides a large image size with excellent contrast and resolution.

OES Pro Nephroscope

WA33036A OP telescope,
30° direction of view,
4 mm channel,
22 Fr. x 230 mm

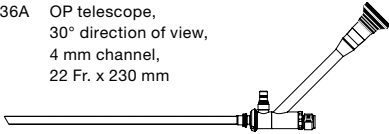


Delivery includes: attachment A3342

WA03200A Light-guide cable,
size S,
plug type,
3 m



A0460
O0332 Light-guide adapter,
for OES 4000 light-guide cables
for Wolf light-guide cables

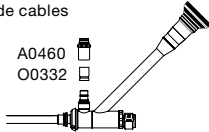

Attachments for Nephroscope

A3342 Attachment,
with automatic valve



WA33039A Attachment,
with stopcock



WA33024A Adapter,
LUS-2 handle
to WA33036A



Use of Fiberscopes

WA33040A Guiding tube,
for cystonephrofiberscope CYF-5/CYF-5A




Continuous Flow Irrigation

Outer sheath,
25 Fr.,
rotatable
WA33037A
WA33035A fixed



WA33025A Obturator, for WA33035A/-37A,
for use of WA33036A
as a compact cystoscope

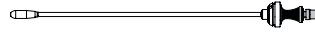


WA33026A Adapter,
OES Pro outer sheath to
Ellik evacuator/bladder syringe



Overflow Irrigation

Irrigation ring,
A22051A 1 stopcock, rotatable
A22052A 1 Luer-lock connector, rotatable
A22053A 2 stopcocks, rotatable
A22054A 1 vertical stopcock, fixed
A22055A 1 vertical Luer-lock connector, fixed

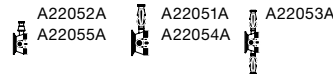

WA33038A Sheath, acc. to Amplatz,
reusable, 28 Fr.



Lumbotome according to Korth

A3328 Knife,
guided, double



Bougies

WA33027A Bougie tubes, set 9–25 Fr.




WA33028A Guiding tube,
for second guide wire,
for WA33027A




Hand Instruments


WA63810A Grasping forceps,
toothed,
3.25 x 400 mm,
with lumen

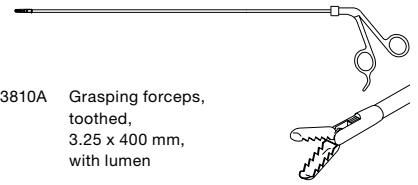


WA63813A Grasping forceps,
3.25 x 400 mm,
with lumen

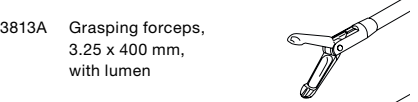


WA63815A Grasping forceps,
3.25 x 400 mm,
toothed

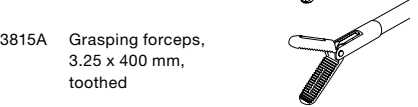

WA63817A Grasping forceps,
three-nail,
3.5 x 400 mm,
with channel 1.1 mm

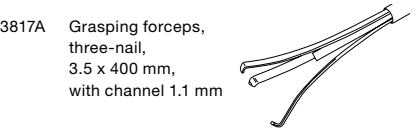


WA63819A Grasping forceps,
three-nail,
3.25 x 400 mm

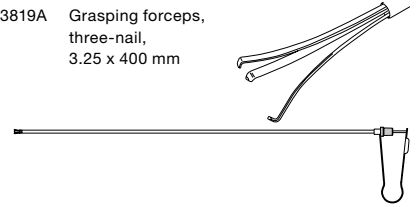


Monopolar Biopsy Forceps

A3334 Biopsy forceps,
4 x 400 mm,
monopolar,
CAM-joint

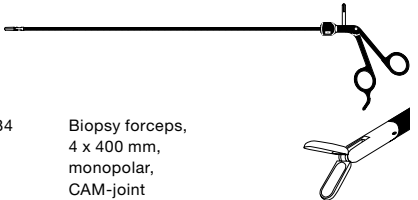

A0358 HF cable,
monopolar,
3.5 m length,
for UES-30/-40,
Erbe ICC International,
and new generation Valleylab HF units
A0355 for UES-10/-20 and old
generation Valleylab HF units,
A0357 for Martin and Berchtold HF units

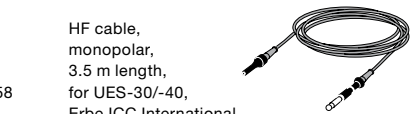

OLYMPUS

Manufactured by
OLYMPUS WINTER & IBE GMBH

Distributed by
OLYMPUS MEDICAL SYSTEMS CORPORATION
OLYMPUS MEDICAL SYSTEMS EUROPA GMBH
OLYMPUS SURGICAL & INDUSTRIAL AMERICA INC.
OLYMPUS AMERICA INC.
KEYMED (MEDICAL & INDUSTRIAL EQUIPMENT) LTD.
OLYMPUS SINGAPORE PTE LTD.
OLYMPUS MOSCOW LIMITED LIABILITY COMPANY
OLYMPUS AUSTRALIA PTY. LTD.
OLYMPUS LATIN AMERICA INC.
OLYMPUS KOREA CO., LTD.
OLYMPUS (BEIJING) SALES & SERVICE, CO., LTD.

Kuehnstraße 61, 22045 Hamburg, Germany

2951 Ishikawa-cho, Hachioji-shi, Tokyo 192-8507, Japan
Wendenstraße 14-18, 20097 Hamburg, Germany
One Corporate Drive, Orangeburg, New York 10962, U.S.A.
Two Corporate Center Drive, Melville, N.Y. 11747-3157, U.S.A.
KeyMed House, Stock Road, Southend-on-Sea, Essex SS2 5QH, United Kingdom
491B, River Valley Road #12-01/04, Valley Point Office Tower, Singapore 248373, Singapore
117071, Moscow, Malaya Kaluzhskaya 19, bld. 1, fl.2, Russia
31 Gilby Road, Mount Waverley, Victoria 3149, Australia
6100 Blue Lagoon Drive, Suite 390 Miami, FL 33126-2087, U.S.A.
8F, Hyundai Marines Bldg., 646-1, Yeoksam-Dong, Kangnam-Gu, Seoul 135-080, Korea
13F - 14F, East Tower, Gongyuan No. 6 Royal Palace, No. 6 Gongyuan Xijie,
Jiangoumennei, Dongcheng District, Beijing 100005, China