

# OLYMPUS<sup>®</sup>

Your Vision, Our Future

Rhino-Laryngo Videoscope

## OLYMPUS ENF TYPE V2

*Slimmer Design for Gentler Insertion—  
Expanding the Possibilities with an Ultra-slim Videoscope*

- High quality imaging

- Smooth insertion with 3.2 mm diameter

- Multi-processor compatibility



3.2 mm ultra-slim

**A bold, bright picture. A slim, sleek design.**  
**Observation and insertion performance have evolved.**

## imaging



### High quality images

The ENF TYPE V2 delivers some of the brightest, clearest images in a videoscope, making observation less strenuous, facilitating a detailed image for diagnosis.




### NBI

Observation with Narrow Band Imaging™ (NBI) is available when used with the EVIS EXERA II™ Universal Platform.



## ultra-slim

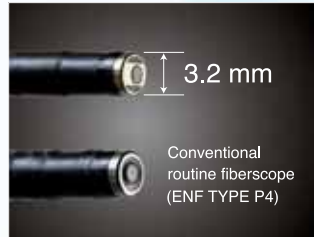


### 3.2 mm diameter

Even thinner than conventional routine fiberscopes, this scope's slim design helps to facilitate insertion while enabling clear observation.




Allows smooth insertion into narrow nasal passages and throat.

## performance

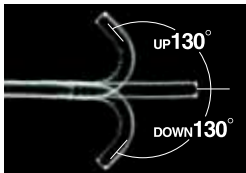

### Flexibility

130°UP/DOWN flexible angulation facilitates less invasive insertion and effective orientation capability.

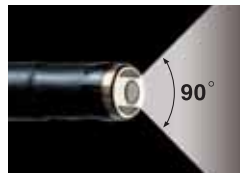

### A Wide FOV

A 90° wide field of view allows a large area to be observed without unnecessary angulation.

### Specifications



<b>Optical system</b>	
Field of View	: 90°
Directions of view	: 0°(Forward viewing)
Depth of field	: 5-50 mm
<b>Insertion tube</b>	
Distal end outer diameter	: 3.2 mm
Flexible tube outer diameter	: 3.4 mm
Working length	: 300 mm
<b>Bending section</b>	
Angulation Range	: Up 130°, Down 130°
<b>Total length</b>	
	: 528 mm
<b>NBI</b>	
	: Available

### Compatible Video Systems



EVIS EXERA II™  
 Universal Platform



VISERA Video System



OTV-SI Compact Office  
 Video System

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

**OLYMPUS®**

OLYMPUS MEDICAL SYSTEMS CORP.  
 Shinjuku Monolith, 3-1 Nishi-Shinjuku 2-chome, Shinjuku-ku, Tokyo 163-0914, Japan  
 OLYMPUS MEDICAL SYSTEMS EUROPA GMBH  
 Postfach 10 49 08, 20034 Hamburg / Wendenstrasse 14-18, 20097 Hamburg, Germany  
 OLYMPUS AMERICA INC.  
 3500 Corporate Parkway, P.O. Box 610, Center Valley, PA 18034-0610, U.S.A.  
 OLYMPUS SURGICAL & INDUSTRIAL AMERICA INC.  
 One Corporate Drive, Orangeburg, New York 10962, U.S.A.  
 OLYMPUS LATIN AMERICA, INC.  
 5301 Blue Lagoon Drive, Suite 290 Miami, Florida 33126-2097, U.S.A.  
 KEYMED LTD.  
 Keymed House, Stock Road, Southend-on-Sea, Essex SS2 5QH, England

OLYMPUS SINGAPORE PTE LTD.  
 491B River Valley Road #12-01/04, Valley Point Office Tower, Singapore 248373  
 OLYMPUS HONG KONG AND CHINA LIMITED  
 Room 1520-1527 Ocean Centre, 5 Canton Road, Tsimshatsui, Kowloon, Hong Kong  
 OLYMPUS (BEIJING) SALES & SERVICE CO., LTD.  
 Room.1202, NCI Tower, A12 Jianguomenwai Avenue Chaoyang District Beijing, 100022, China  
 OLYMPUS MOSCOW LIMITED LIABILITY COMPANY  
 117071, Moscow, Malaya Kaluzhskaya 19, bld. 1, 1L2, Russia  
 OLYMPUS AUSTRALIA PTY LTD  
 31 Gilby Road, Mount Waverley, VIC., 3149, Australia

www.olympus.com